Istruzioni per l'uso
Instructions for use
Instrucciones para el uso

SEE-TEST

Articolo
Article
Artículo:
260/S

ITALIANO
ENGLISH
ESPAÑOL

ML

1 - 12

ITALIANO

14 – 23

ENGLISH

25 – 34

ESPAÑOL

36 - 45

1.	Avvertenze generali e informazioni al destinatario General warning and information for the user Advertencias generales e informaciones para el destinatario
1.1	Dichiarazione di Conformità / Declaration of Conformity / Declaración de Conformidad

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GUARANTEE CONDITIONS

- 1) duration: the product is covered by a 12-month (twelve) guarantee.
- 2) exclusion: the guarantee excludes:
 - a) Faults originating from non-performance of scheduled maintenance and due to negligence or improper use of the product by the user;
 - b) Periodic checks and maintenance;
 - c) Repair or replacement of parts subject to wear that are fragile or have an unpredictable life, unless they prove faulty at the time of their delivery;
 - d) Faults caused by personnel or originating from technical personnel transfer and transportation;
 - e) Faults or damage originating from improper use or operating errors;
 - f) Faults or damage originating from contamination in the water and air infeed systems, extraordinary chemical agents or electrical events;
 - g) Faults or damage originating from the use of detergents, disinfecting agents, sterilizing fluids or processes that are not described explicitly in the instruction manual for operation and maintenance;
 - h) Natural color change of plastic components.
- 3) installation and testing report: an essential prerequisite for guarantee validity is the return of the installation and testing report after all its sections have been filled in by the user and the fitter. The report must be returned within 15 days after installation, otherwise the guarantee becomes null and void.
- 4) limits: the guarantee entitles its holder to repairs or replacement of faulty components free of charge. The right to replace the entire machine is excluded. As for components either applied or added to tecno-gaz products or by third parties, having their own guarantee certificate, the conditions, limits and exclusions indicated in the relevant certificates apply.

Under no circumstances shall carrying out of one or more than one repair during the guarantee term change the guarantee expiry date.
- 5) claims: in the event of the purchaser making a claim relating to the application of the guarantee and the quality or conditions of the delivered equipment, the said purchaser shall in no way suspend and/or delay payments.
- 6) disputes: if a dispute arises on the application and interpretation of this certificate of guarantee, the court of parma (italy) shall be the court of jurisdiction, regardless of the location in which the purchase agreement was stipulated.
- 7) derogations: any special derogations from the guarantee terms described in this section shall not imply the granting to the purchaser of any right and shall be considered as granted for a specific case only.
- 8) miscellaneous: for any other issue, which is not covered by this certificate of guarantee, the rules described in the civil code and existing relevant laws in the republic of italy apply.

LOSS OF THE GUARANTEE

The guarantee becomes null and void in the event that:

- a) the equipment presents damage caused by a fall, exposure to flames, liquid spillage, natural events, bad weather conditions or by any other cause that is not due to manufacturing defects;
- b) installation has not been performed in line with TECNO-GAZ instructions and by authorized personnel;
- c) the equipment is repaired, modified or tampered with by the Purchaser or by non-authorized third parties;
- d) when asking for an intervention under guarantee, the product serial number has been removed, deleted, counterfeited, etc.;
- e) the installation and testing report is not returned suitably filled in and signed within 15 days from installation;
- f) the Purchaser suspends and/or delays payment of amounts due, for any reason whatsoever, relating to the purchase and/or maintenance of the equipment;
- g) scheduled periodic maintenance or any other instruction envisaged in the manual for operation and maintenance is not carried out.

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DELIVERY TERMS

1. All goods delivered must be packaged in their original packaging. If the goods are delivered in a non-genuine packaging, a charge will be applied to restore the original packaging.
2. All goods must be delivered CARRIAGE PAID.
3. The client shall ask a written authorization to TECNO-GAZ sales dept. before returning the goods by means of the "Form for Client's Returned Goods Management" (CM-M-002) in compliance with the requirements of TECNO-GAZ Quality Management System (ref. CM-P-003, CM-I-009).
4. The document accompanying the returned goods shall show:
 - a. the TECNO-GAZ product ode/item no.;
 - b. the number of the TECNO-GAZ purchase invoice or shipment document;
 - c. the product serial number and batch number;
 - d. a description of the goods returned;
 - e. the number of the Form for Client's Returned Goods Management" (CM-M-002) or a copy of the same attached.
5. Returned goods are considered as accepted only after they have been examined. If they prove non-conforming, they may be sent back to the consignor at TECNO-GAZ full discretion.
6. Charges will be applied to restore and test the returned product.
7. Parts replaced under guarantee must be returned, carriage paid, to TECNO-GAZ.
8. If the parts are not returned, their cost will automatically be debited to the recipient.
9. TECNO-GAZ does not accept goods returns from end users.
10. Goods sent to TECNO-GAZ for repair are managed as returned goods.

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E-Mail : info@tecnogaz.com

<http://www.tecnogaz.com>

For technical service please apply to the dealer that invoiced the device, or directly to:

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
Par 09

Schemistica

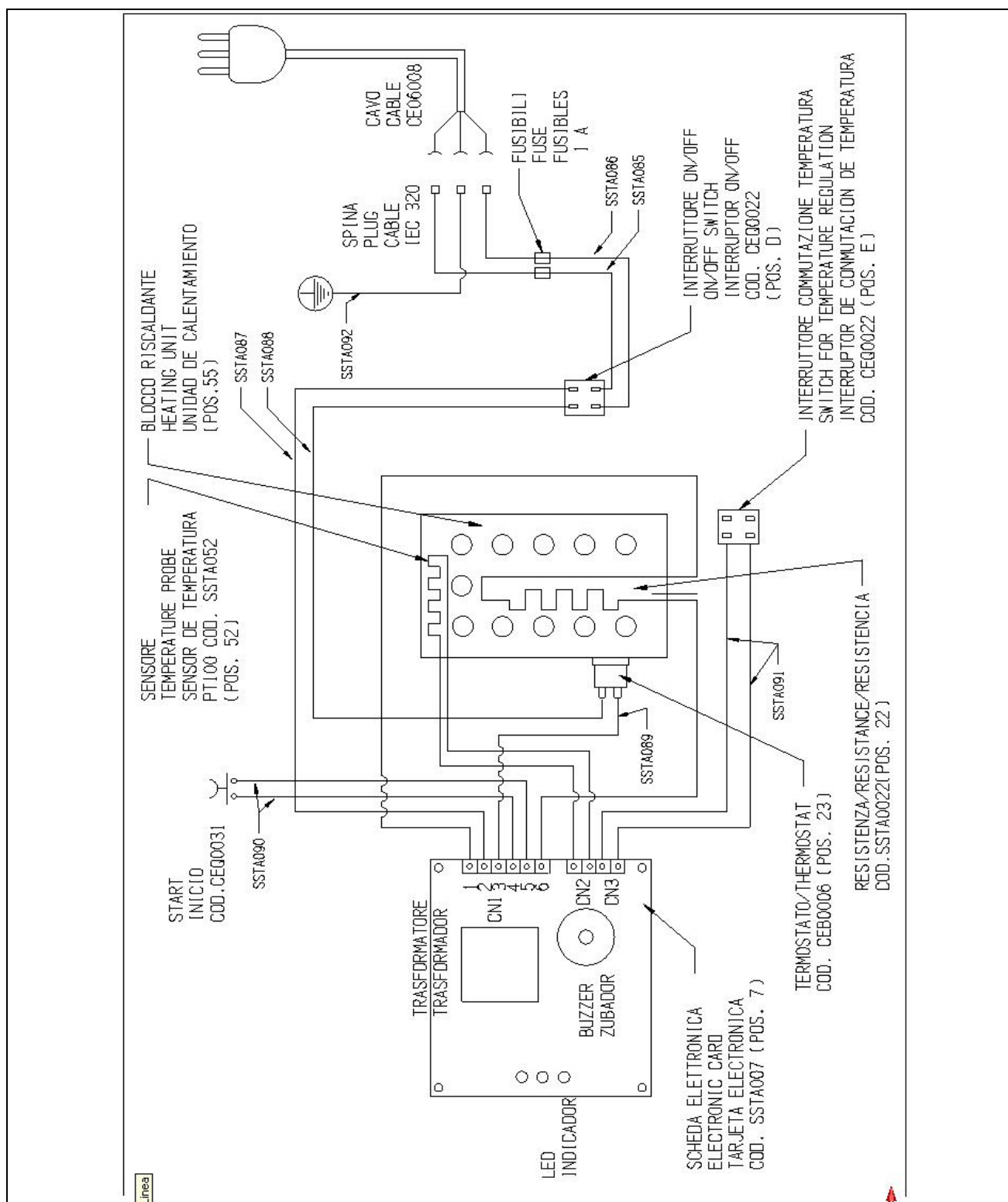INDEX		
PAR	DESCRIPTION	PAG
01	INFORMATION TO THE USER	1-12
02	INTRODUCTION PRODUCT	25-26-27
03	TECHNICAL FEATURES	28
04	PACKING – TRANSPORT - STORAGE 4.1 INSTALLATION 4.2 DEMOLITION	29 30 30
05	WORKING AND USE 5.1 WORKING 5.2 RELIABILITY OF THE TEST VALIDITY	31 31 31
06	INSTRUCTION FOR USER 6.1 WORKING FOR STEAM AT 57°C 6.2 WORKING FOR GAS AT 37°C	32 32 32
07	ORDINARY AND EXTRAORDINARY MAINTENANCE	33
08	ANALYSIS AND RESEARCH OF FAILURES	33
09	DIAGRAM	34

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The See-Test (PICT.1) is a device which was developed for the incubation of biological indicators at single temperature.


The incubator is basically made up of a cabinet/support, a thermoplastic receptacle, an aluminium heating unit (with 15 seats for test-tubes/biological indicators) and one electronic card for the automatic regulation of the heating unit temperature.

On the back of the device you can find:

- Main switch ON/OFF (POS.D-FIG.2A)
- Temperature regulation switch (POS.E – FIG.2A)
- Networking plug (POS.F – PICT.2A)

ATTENTION: DO NOT TOUCH THE SURFACE OF THE HEATING UNIT WITH WORKING MACHINE .

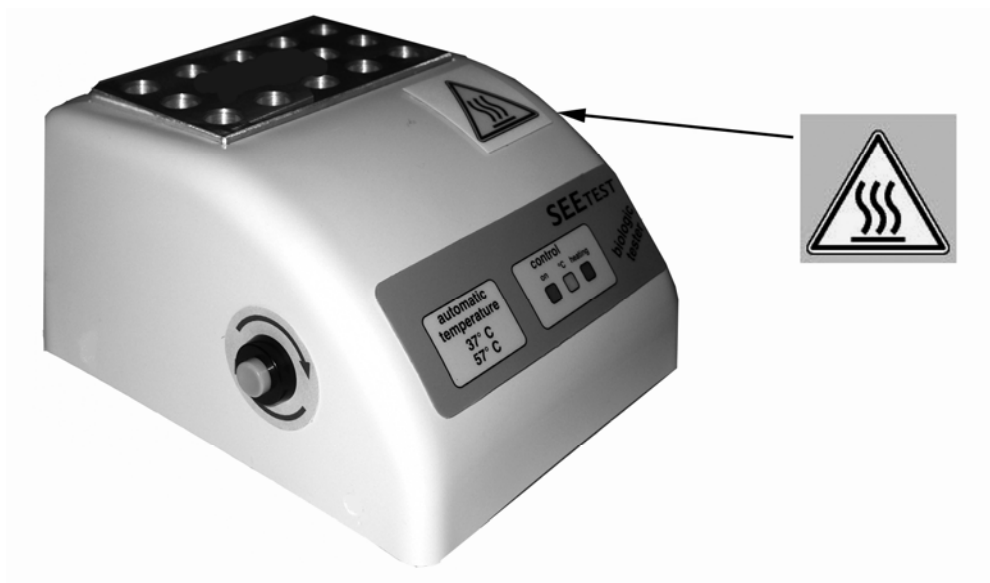
PICT. 1


PICTOGRAM

ATTENTION:


**HOT
SURFACE**

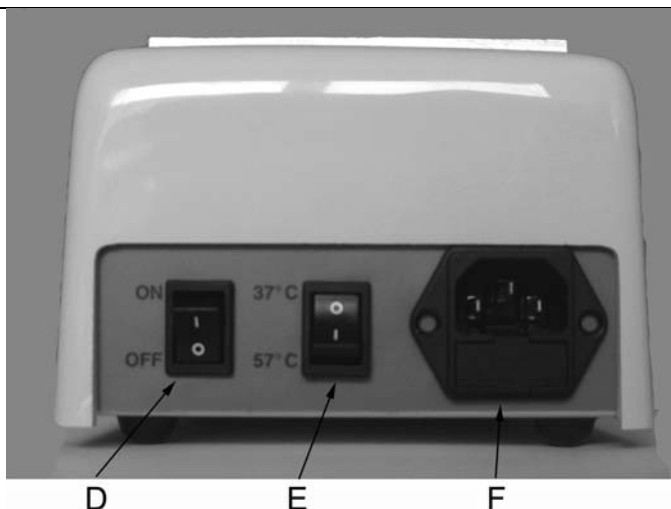PICT. 2A

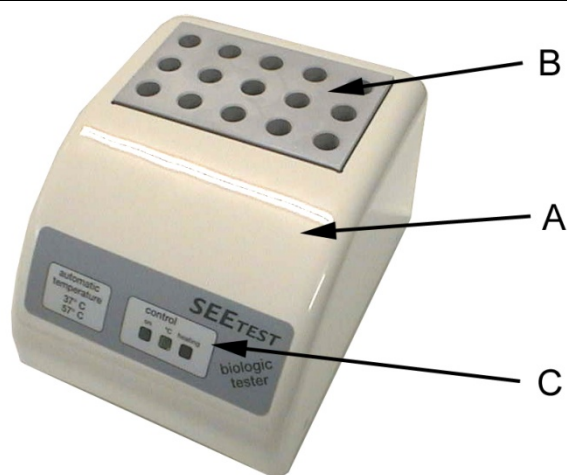PICT. 2B

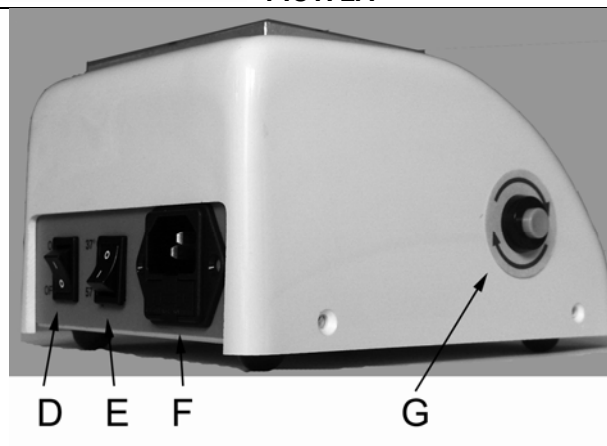PICT. 2C

PICT	POS	DESCRIPTION
2B	A	THERMOPLASTIC COVER
2B	B	HEATING UNIT WITH N°14 SEATS FOR TEST-TUBES
2B	C	LUMINOUS INDICATOR OF TEMPERATURE
2A-2C	D	MAIN ON/OFF SWITCH
2A-2C	E	SWITCH FOR TEMPERATURE REGULATION
2A-2C	F	FUSE BLOCK PLUG
2C	G	START BUTTON

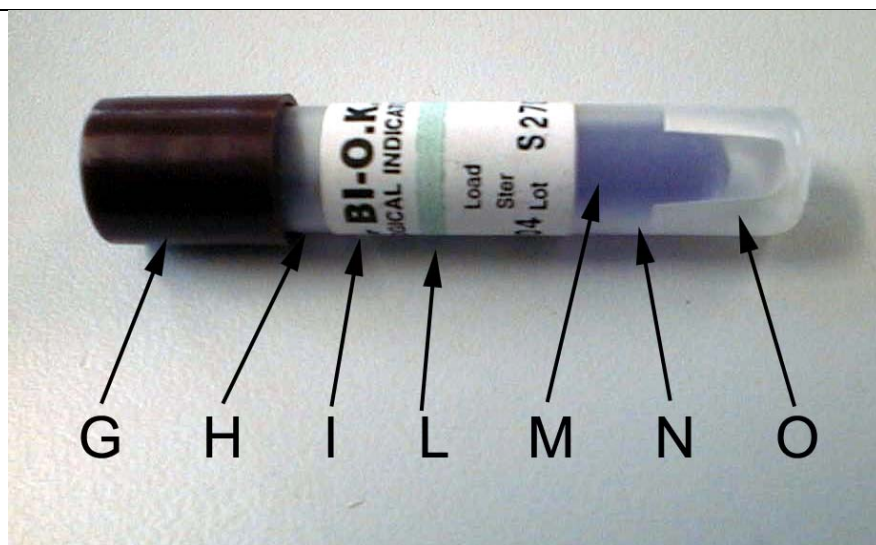
FIG 3

PICT.	POS	DESCRIZIONE
3	G	LID
3	H	FILTER
3	I	LABEL
3	L	TEST-TUBE
3	M	GLASS BULB
3	N	CULTURE MEDIUM
3	O	PAPER STRIPE


FIG 4

Par 03	Technical features
---------------	---------------------------

MECHANICAL FEATURES OF THE DEVICE	
Working room temperature	17°C ÷ 27°C
Sizes (X,Y,Z)	140 x 100 x 170 mm
Weight	1,8 kg
Body material	Fe37 steel RAL 9002 painted
Cover material	Shrunked ABS
Heating unit material	Anodised aluminium
Packaging sizes (X,Y,Z)	205 x 160 x 305 mm
ELECTRICAL FEATURES OF THE DEVICE	
Supply voltage	230 V~ ± 10%
Number of phases	1
Supply frequency	50/60 Hz
Max. absorbed power	50 W
Protection degree	PE
Insulation class	1
Length of the supply cable l= 2200 mm	2 x 1 sqm + GROUND
Fuses	1A 250V
Electrical connector (plug)	ITALY PLUG 10 A + EARTH CONNECTION IEC 320

X= Width Y=Height Z= Depth

Par 04	Packing – transport – storage – installation – demolition - disposal
--------	--

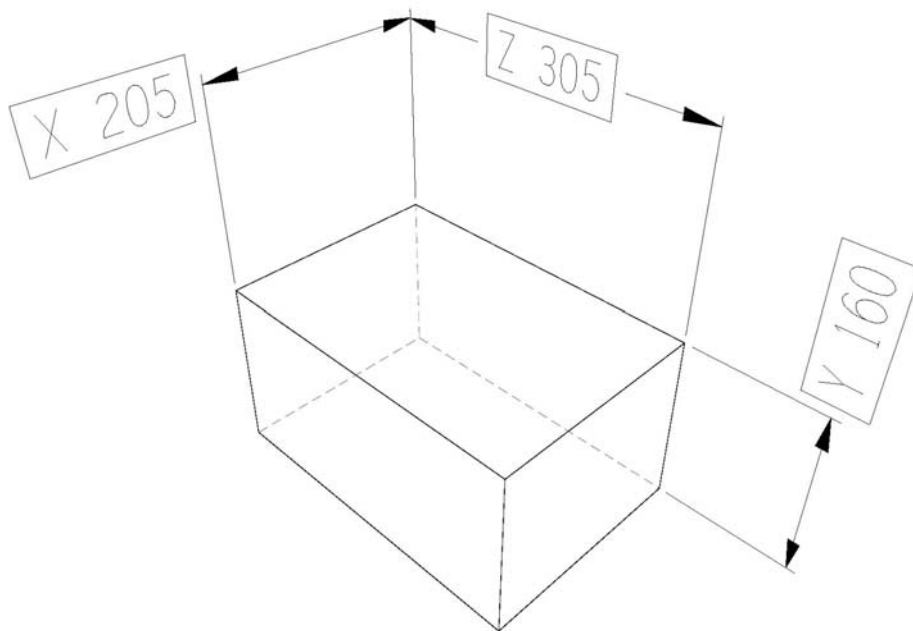WSTM001
The specified sizes
are expressed in
millimetres (mm)

SEE-TEST is packed as follows: inside a cellophane bag, wrapped in a foam polystyrene cover, placed within a corrugated cardboard box.

Non-sterile package.

The package and the device are delicate, handle with care, make sure they are not bumped and shaken during transport operations, do not turn upside-down.

The packed equipment must be kept in a dry place and at a temperature between -10° and +40°.

It is allowed to place maximum 3 machines on top of one another.

THE PACKAGING MATERIALS MUST BE PRESERVED DURING THE WHOLE WARRANTY PERIOD.
TECNO-GAZ WON'T ACCEPT RETURNS WITHOUT ORIGINAL PACKAGING MATERIAL.

4.1	Installation
------------	---------------------

Installation is a fundamental operation for the machine's proper use and working.
Here are the operations to follow for a correct installation:

Do not keep the SEE-TEST biological incubator in rooms with potential risk of explosion.

Remove the machine from its packaging and check its condition.
Notify TECNO-GAZ in case the machine bears signs of damage caused by transport.

The packaging materials must be preserved during the whole warranty period.
TECNO-GAZ won't accept returns without original packaging.

Locate the machine on a flat surface, away from sources of heat (e.g. autoclaves, bunsen, pre-heating furnaces, etc.) at a maximum range of 1,5 metres from the outlets. These outlets must be handy.

Adapt the height of the work surface to the operator (we suggest 90 cm).

The led lights must be clearly visible to the user.

Connect the plug to an outlet which complies to the features described in the technical specifications.

4.2	Scrapping
------------	------------------

SEE-TEST is manufactured using plastic materials (polystyrene and ABS), ferrous materials and electronic components.

For scrapping and disposal refer to the laws in force in the country where the machine is used and then contact the authorized firms.

DO NOT DUMP THE MACHINE AMONG URBAN RUBBISH

Par 05	Working and application
---------------	--------------------------------

SEE-TEST is a device for medical and/or dental use.

TECNO-GAZ will not be held responsible for damage to the machine or harm to people caused by the improper use of the machine or by the user's unskilfulness.

ROOM TEMPERATURE +10 / +30 °C

5.1	Working
------------	----------------

1. Select the temperature cycle on the temperature selector switch (POS. E – PICT. 2A)

Position 0 corresponds to temperature cycle 37°C

Position 1 corresponds to temperature cycle 57°C

2. Select position 1 on the general switch (POS. D – PICT. 2A), push START button (POS. G – PICT. 2C): "ON" and "HEATER" led lighten. "HEATER" led keep lighten till reaching the selected temperature, then "°C" led lighten and a sound signal advise of the reaching of the selected temperature. "°C" led also means that the selected temperature keep even. When the temperature decrease, due to the natural cooling of the heating element, the resistance begins to work. "ON" led always lighten, this means that the device is working.

5.2	Reliability of the test validity
------------	---

As already said, "ON" led means the device regular working, in case of lack of electrical energy, or broken control system, the device switch off permanently and led keep turned off, also after the reset of the electrical energy. In this way the user, incubation time past, can verify, if led are switched on, that the test is valid. Otherwise, led switch off means an electrical energy interruption or a broken control system. The test is not valid.

Par 06	Instruction for user
6.1	Working for steam at 57°C

1. Identify the biological indicator to be treated.
2. Insert the biological indicator into the sterilizer.
3. Perform the sterilization cycle.
4. Remove the biological indicator from the sterilizer.
5. Let the biological indicator cool for at least 10 minutes.
6. Check the colour change to be sure that the sterilization has been properly carried out because the chemical indicator changes colour when it is subjected to steam.
7. Break the glass bulb of the biological indicator which is inside the thermoplastic container by exerting pressure with the supplied pliers, until the bulb breaks. Perform this operation just a minute before inserting the indicator in the incubator. The indicator must be inserted in the incubator within and not beyond 1 hour from the moment it was removed from the sterilizer.
8. Insert the biological indicator in the incubator's heating unit.
9. Check the biological indicator at regular intervals to control the changing of colour (every 18 hours, 24 hours, 36 hours).
10. Incubate the biological indicator for at least 48 hours using as a reference an indicator not treated in a sterilisation cycle.
11. For the detailed procedure carefully read the instructions for use which are inside the packaging of the B.I.O.K. biological indicators (the biological indicators can be supplied only on demand).

WARNING: "SEE-TEST" is projected to incubate the B.I.O.K. biological indicators (supplied by Propper). TECNO-GAZ will not be held responsible for the results obtained with other indicators.

6.2	Working for gas at 37°C
------------	--------------------------------

1. The process is similar to that for the steam working.
2. Stick rigorously to the technical specifications of the biological indicator which is used.

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The machine does not need any particular maintenance apart from the periodic cleaning of the heating unit which is necessary to avoid the deposit of substances inside the seat of the biological indicators. Use a wet cloth for cleaning. Do not use aggressive cleansers.


DO NOT DO MAINTENANCE OR INSIDE REPAIR TO THE DEVICE.
 DISMANTLE THE DEVICE MEANS BREAKING IT

Par 08	Analysis and research of failures
---------------	--

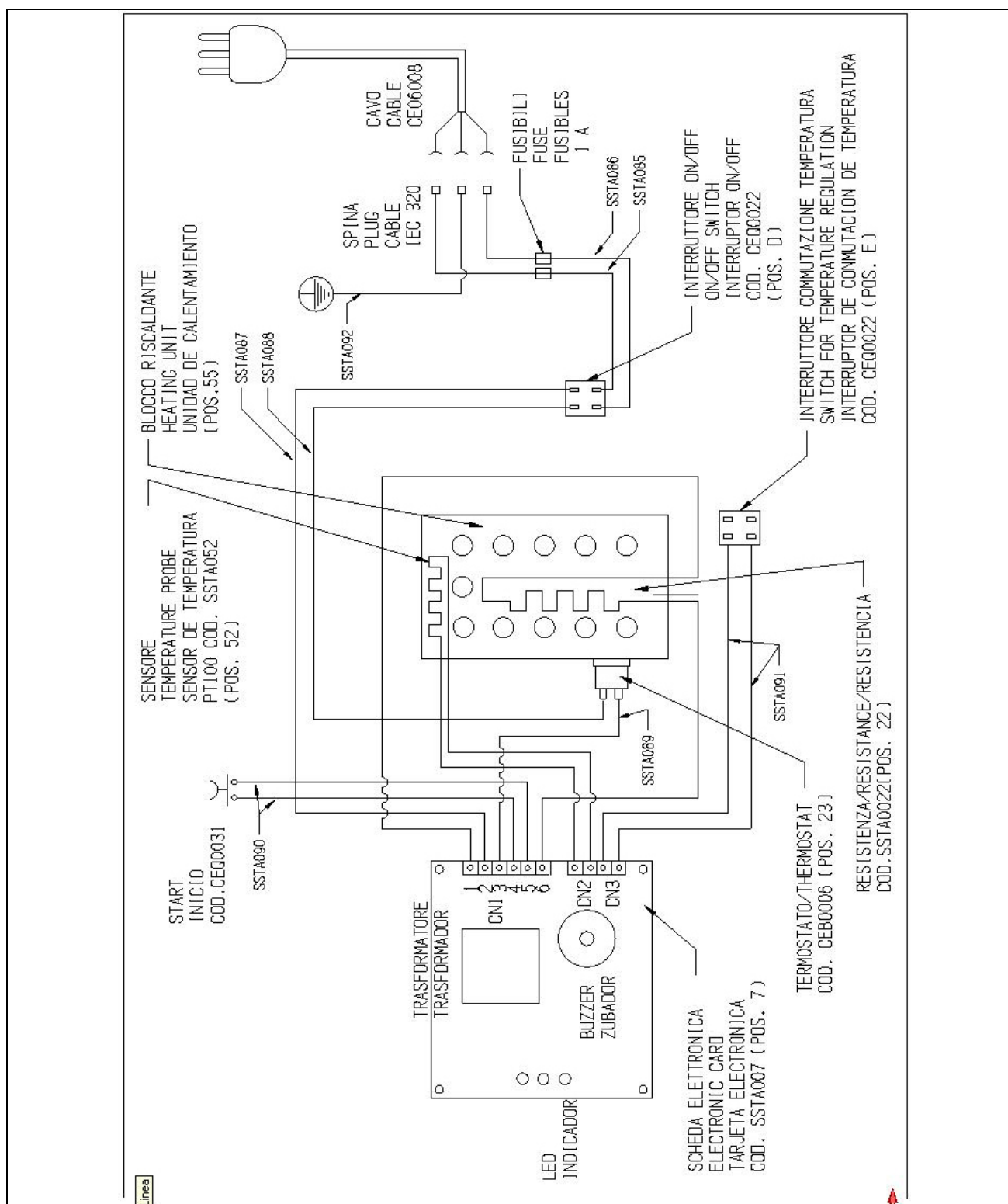N°	Difetto	Causa	Intervento
1	Pushing start button see-test does not work	There is no power supply	Control plug fit in correctly
2	With led switch on see-test does not reach the selected temperature	Anomalies in the electronic card or in the thermostat	Contact Tecno-gaz s.p.a. for technical assistance
3	Led switch off during the test	Electrical energy interruption	Push start button, if led don't switch on, control cable and power plug. If led switch on, remove tubes and keep device switch on for about 30 minutes. If device don't switch off automatically, means that the previous stopping was due to a temporary lack of electrical energy. Otherwise if device switch off, means that the heating unit work over the selected temperature and the safety device become effective. In this case contact Tecno-gaz s.p.a.
		Inside failure	Contact Tecno-gaz s.p.a. for technical assistance

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TECNO-GAZ S.p.A.